


# SafeSheath® Ultra Lite

HEMOSTATIC VALVE PEEL AWAY INTRODUCER SYSTEM

# LATEX-FREE

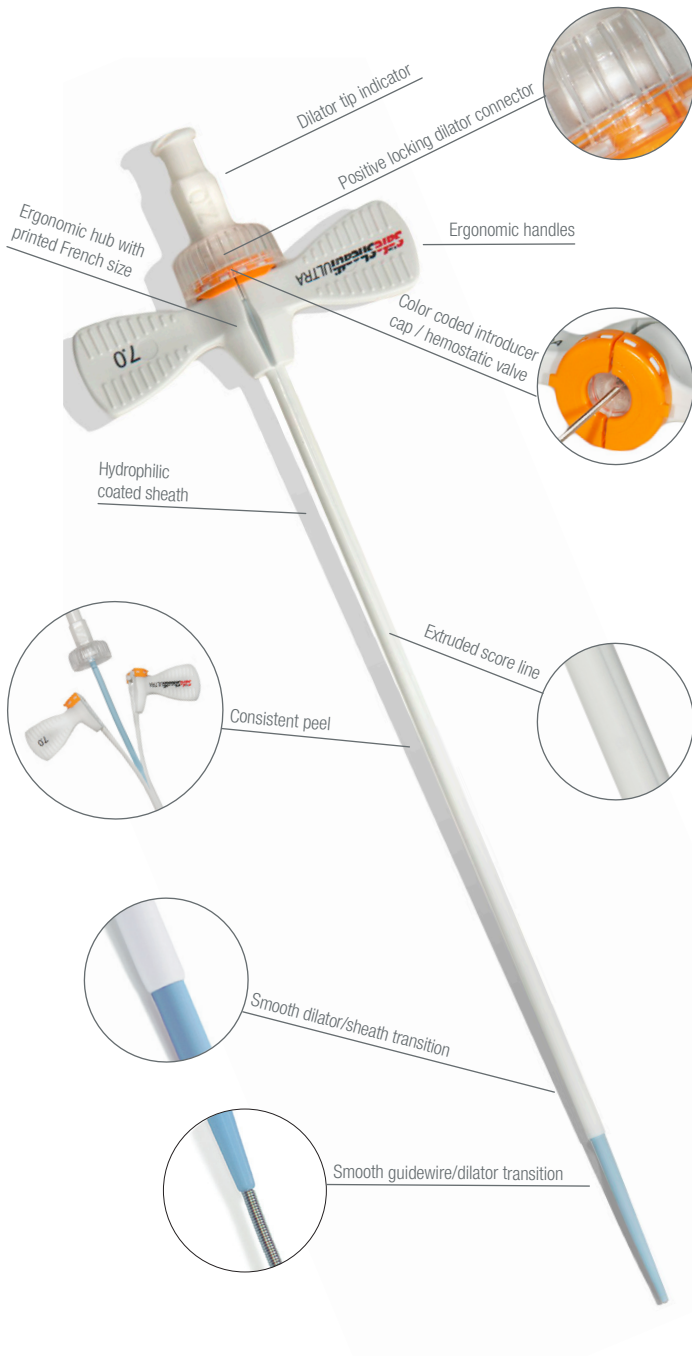

Port and central venous catheters require an introducer with easy insertion properties. Before now, that meant the possibility of significant bleeding or air embolism.

Introducing the **SafeSheath® Ultra Lite Hemostatic Valve Peel Away Introducer System**. Incorporating a specially engineered valve membrane, the **Ultra Lite** provides reduced insertion forces during procedures.

## SafeSheath® Ultra Lite Hemostatic Valve Peel Away Introducer System Features:

- Specially designed membrane allows effortless insertion and removal of ports and central venous catheters
- Extruded score line peel-away sheath for clean, even peeling during removal
- Ergonomically designed, easy splitting breakaway hub
- French size printed on handle and indicated by color-coded introducer cap
- Positive locking dilator connector
- Decreased risk of blood loss and air embolism
- Standard 13 cm length with a wide range of French sizes available

Model / Catalog #	French Size	Sheath Length	Dilator Length
ASUP051318	5 F	13 cm	18 cm
ASUP061318	6 F	13 cm	18 cm
ASUP071318	7 F	13 cm	18 cm
ASUP081318	8 F	13 cm	18 cm
ASUP851318	8.5 F	13 cm	18 cm
ASUP091318	9 F	13 cm	18 cm
ASUP101318	10 F	13 cm	18 cm
ASUP111318	11 F	13 cm	18 cm

All SafeSheath Ultra Lite Introducer Systems include sheath and dilator. The H. Lee SafeSheath® hemostatic splittable sheath is protected by Patent Insurance through 2011. USA Patents 5,125,904 and 5,312,355. Canadian Patent 2,112,394. European Patent EP0593685.

\* SafeSheath is a registered trademark of Pressure Products, Inc.

Introducers supplied sterile. For single use only. Available bulk upon request.

Standard lengths are noted above. Other lengths available upon request.

For additional information e-mail us at [AdelanteSafeSheathUltraLite@oscor.com](mailto:AdelanteSafeSheathUltraLite@oscor.com)

## SafeSheath® Ultra Lite Sterile Kit Contains:

- One SafeSheath® Ultra Lite peel-away sheath
- One dilator
- One guidewire with a flexible straight and J tip with straightener
- One needle with protector
- One syringe



Available in other kit configurations

