

# Petite™

THE WORLD'S THINNEST PASSIVE  
FIXATION BIPOLAR PACING LEAD


# LATEX-FREE

## Petite™ Series Cardiac Pacing Lead

The Petite™ lead represents Oscor's next generation of permanent implantable cardiac pacing leads. The revolutionary design of the Petite™ lead enables a minimal lead body diameter while maintaining maximum reliability. Available in both straight and pre-formed J configurations, Petite™ leads are ideal for either ventricular or atrial placement. With a lead body diameter of 4.8 F and the ability to pass through a 7 F introducer, the Petite™ lead is one of the thinnest bipolar endocardial pacing leads available today. The uniquely designed fixation tines allow for safe and easy fixation of the lead and maximize stability in the heart.

### Petite™ Passive Fixation Lead Features:

- Straight and pre-formed J configurations available for ideal placement in the ventricle or atrium
- Hexafilar coil consists of two conductors connected to the anode and two conductors connected to the cathode
- Each set of conductors is separated by a passive, non-electrically active filar for additional safety
- Polyurethane outer insulation is highly reliable for long term implantation while maintaining lead flexibility
- Polaris™ (Iridium Oxide Pt/Ir) coated distal tip featuring low polarization which may lower acute and chronic pacing thresholds
- Available with steroid eluting ring, dexamethasone sodium phosphate, which has been associated with reduced acute and lower chronic thresholds. *The steroid eluting ring is not available for US market.*

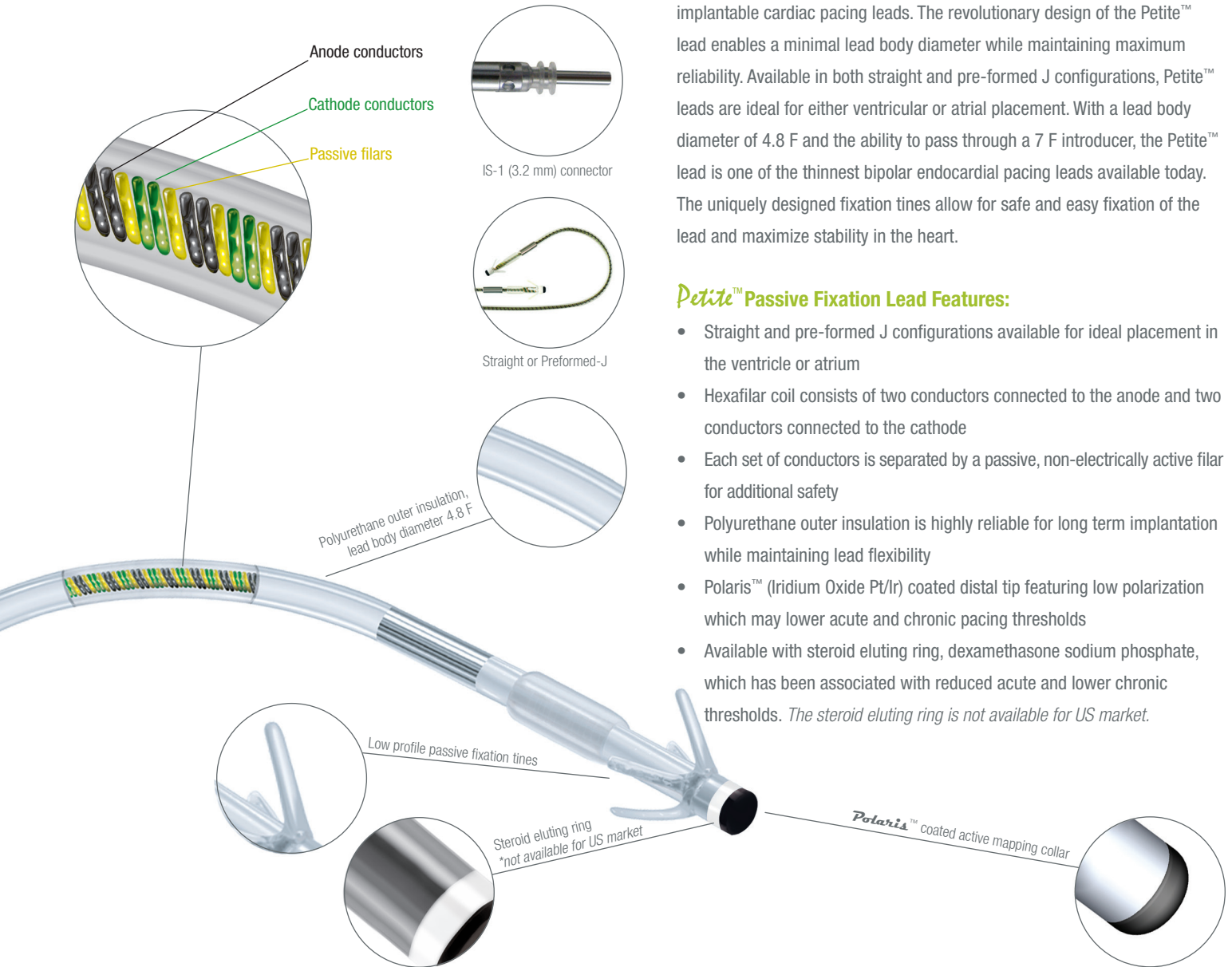

Model No. USA	Model No. non-USA	Placement of Stimulation	Lead Length	Lead Body Diameter	Polarity	Outer Material	Electrode Spacing	Connector Type
Petite 52 RB	Petite 52 ERB	Ventricle	52 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 58 RB	Petite 58 ERB	Ventricle	58 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 70 RB	Petite 70 ERB	Ventricle	70 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 78 RB	Petite 78 ERB	Ventricle	78 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 44 RJB	Petite 44 ERJB	Atrium	44 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 48 RJB	Petite 48 ERJB	Atrium	48 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 52 RJB	Petite 52 ERJB	Atrium	52 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Petite 70 RJB	Petite 70 ERJB	Atrium	70 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)

Leads supplied sterile. For single use only.

Standard lengths are noted above. Other lengths available upon request.

4.8 F Polyurethane Lead body

7 F Introducer