


Physique®

THE WORLD'S THINNEST EXTENDABLE
RETRACTABLE ACTIVE FIXATION LEAD

LATEX-FREE

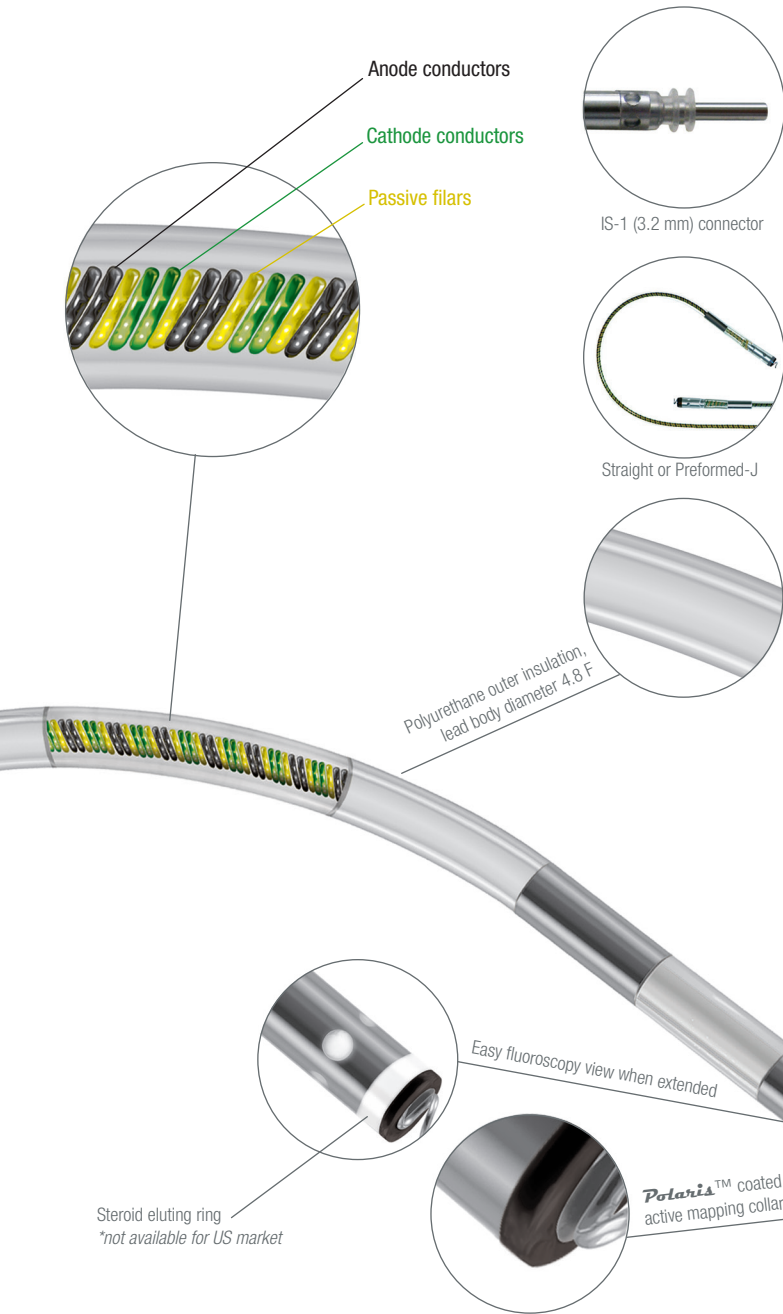

Physique® Series Cardiac Pacing Lead

The Physique® lead represents Oscor's next generation of small profile endocardial permanent pacing leads with extendable/retractable active-fixation helix. With a lead body size of only 4.8 F and a lead tip size of 7 F, the Physique® lead is the smallest extendable/retractable bipolar pacing lead available today.

Physique® Active Fixation Lead Features:

- Straight and pre-formed J configurations available for ideal placement in the ventricle or atrium
- Helix (1.5 mm) allows for easy fixation of the screw in the myocardium
- Screwdriver stylet allows helix to be extended/retracted in only 2-3 turns
- Hexafilar coil consists of two conductors connected to the anode and two conductors connected to the cathode
- Each set of conductors is separated by a passive, non-electrically active filar for additional safety
- Polyurethane outer insulation is highly reliable for long term implantation while maintaining lead flexibility
- The distal portion between the cathode and anode is highly flexible to help prevent perforation of the myocardium
- Equipped with four side-view holes which are clearly visible under fluoroscopy when the screw is fully extended
- Polaris™ (Iridium Oxide Pt/Ir) coated distal tip featuring low polarization which may lower acute and chronic pacing thresholds
- Available with steroid eluting ring, dexamethasone sodium phosphate, which has been associated with reduced acute and lower chronic thresholds. *The steroid eluting ring is not available for US market.*

Model No. USA	Model No. non-USA	Placement of Stimulation	Lead Length	Lead Body Diameter	Polarity	Outer Material	Electrode Spacing	Connector Type
Physique 52 RB	Physique 52 ERB	Ventricle	52 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Physique 58 RB	Physique 58 ERB	Ventricle	58 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Physique 44 RJB	Physique 44 ERJB	Atrium	44 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Physique 48 RJB	Physique 48 ERJB	Atrium	48 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)
Physique 52 RJB	Physique 52 ERJB	Atrium	52 cm	4.8 F	Bipolar	Polyurethane	11 mm	IS-1 (3.2 mm)

Leads supplied sterile. For single use only.

Standard lengths are noted above. Other lengths available upon request.

4.8F Polyurethane Lead body

7F Introducer