

 **oscor**[®]

Destino[™]



Destino[™]

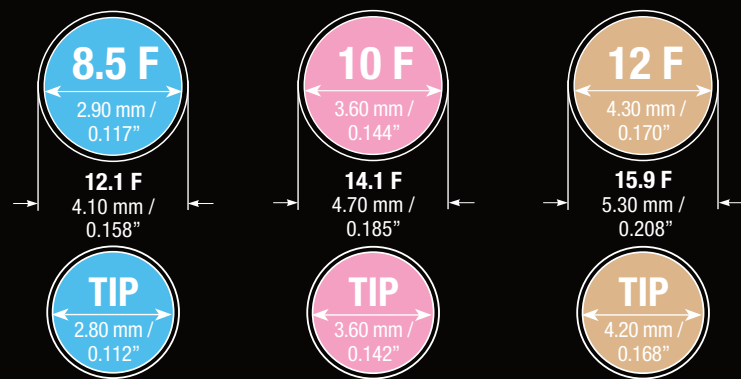
STEERABLE GUIDING SHEATH

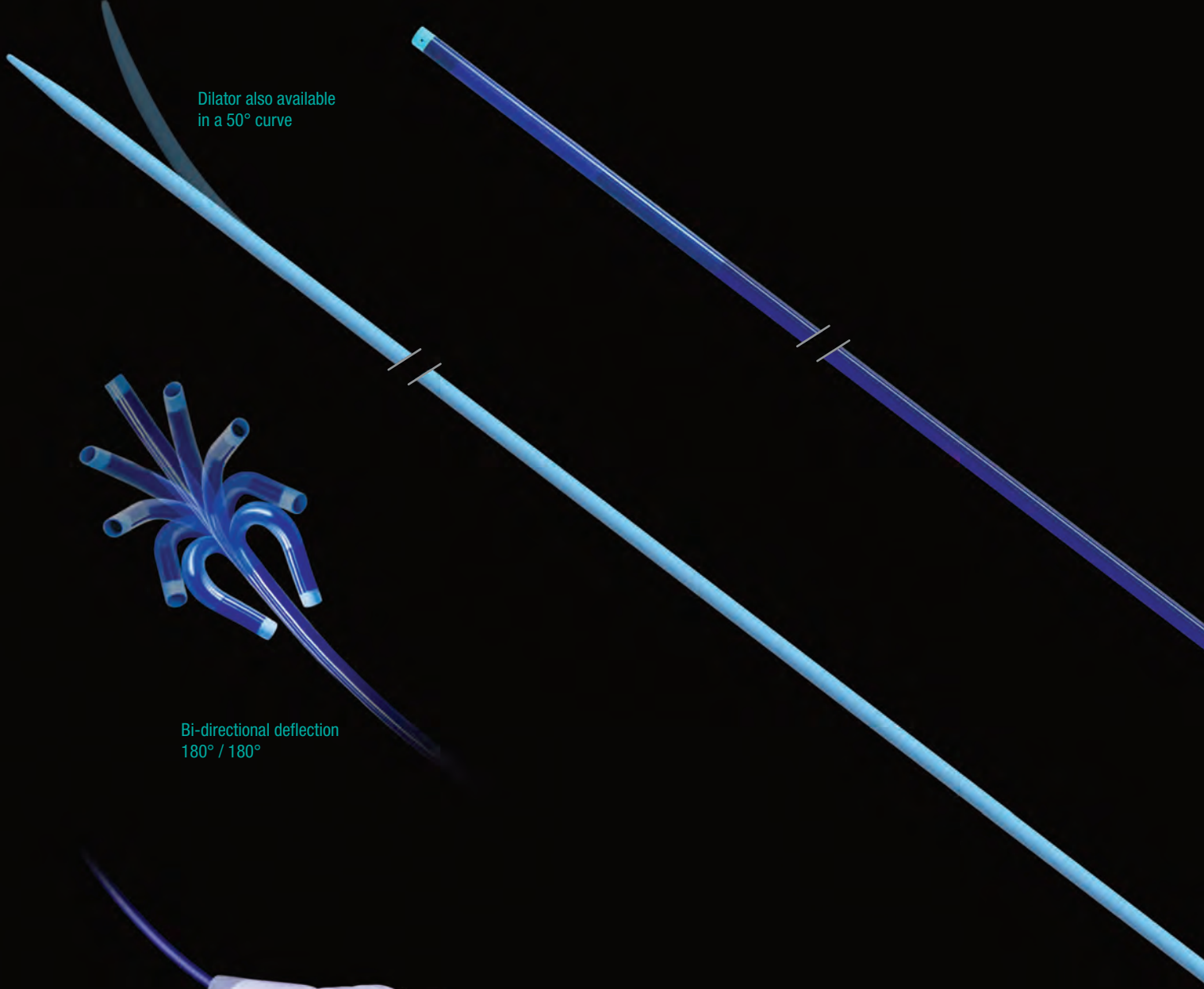


The Oscor *Destino* is a next generation bi-directional steerable guiding sheath for the introduction and placement of diagnostic and therapeutic devices into the human vasculature, including but not limited to intracardiac, transseptal, and other placements.

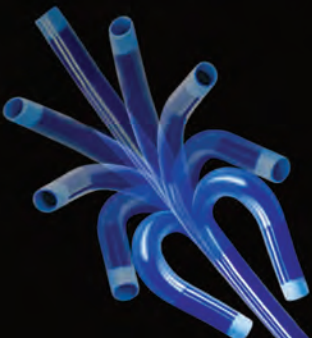
Featuring a next generation hemostatic valve and braided FlexSteer® Shaft Technology, the *Destino* is the only available steerable sheath in a size range of 8.5 F, 10 F, and 12 F with true bi-directional tip deflections of 180°.

Destino Oversized Shaft Inner-lumen and Outer Diameter:






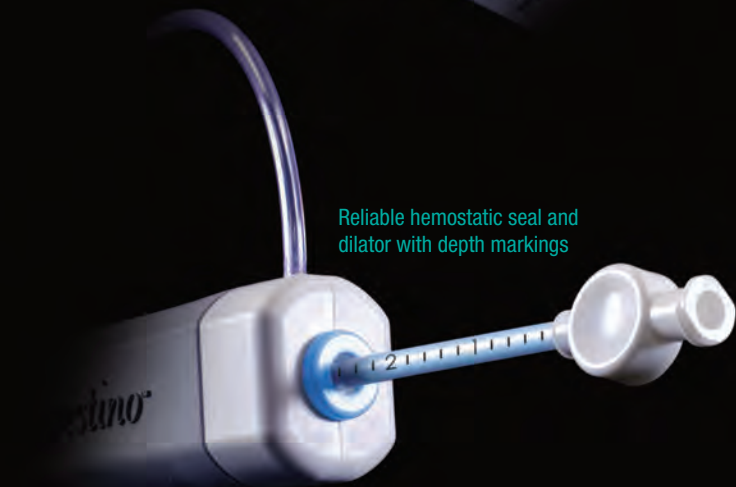
Dilator also available
in a 50° curve



Bi-directional deflection
180° / 180°



Easy locking ergonomic
steerable handle



Reliable hemostatic seal and
dilator with depth markings

The *Destino*® Steerable Guiding Sheath Features:

- Bi-directional steering allows **dual deflection** for maximum maneuverability
- **Ergonomic handle** with curve diameter and French size indicator provides precise control and more comfort
- **Reliable valve** provides effective hemostasis and seals down to a 0.014 inch guidewire
- Maximum compatible guidewire 0.038"/0.97 mm
- Dilator tip construction and depth markings ensures **exact needle exposure**

The *Destino*® Steerable Guiding Sheath Features:

- Soft, atraumatic tip reduces potential for trauma during transseptal punctures and provides smooth transition to catheter
- Radiopaque tip marker enhances fluoro visibility
- Steerable hydrophobic coated sheath improves access to hard-to-reach sites and eliminates need to change out sheaths to reach desired position
- Braided Flexsteer™ shaft provides exceptional torqueability, pushability, and kink resistance



Advanced braided shaft technology



Smooth tip to dilator transition / two sheath side holes

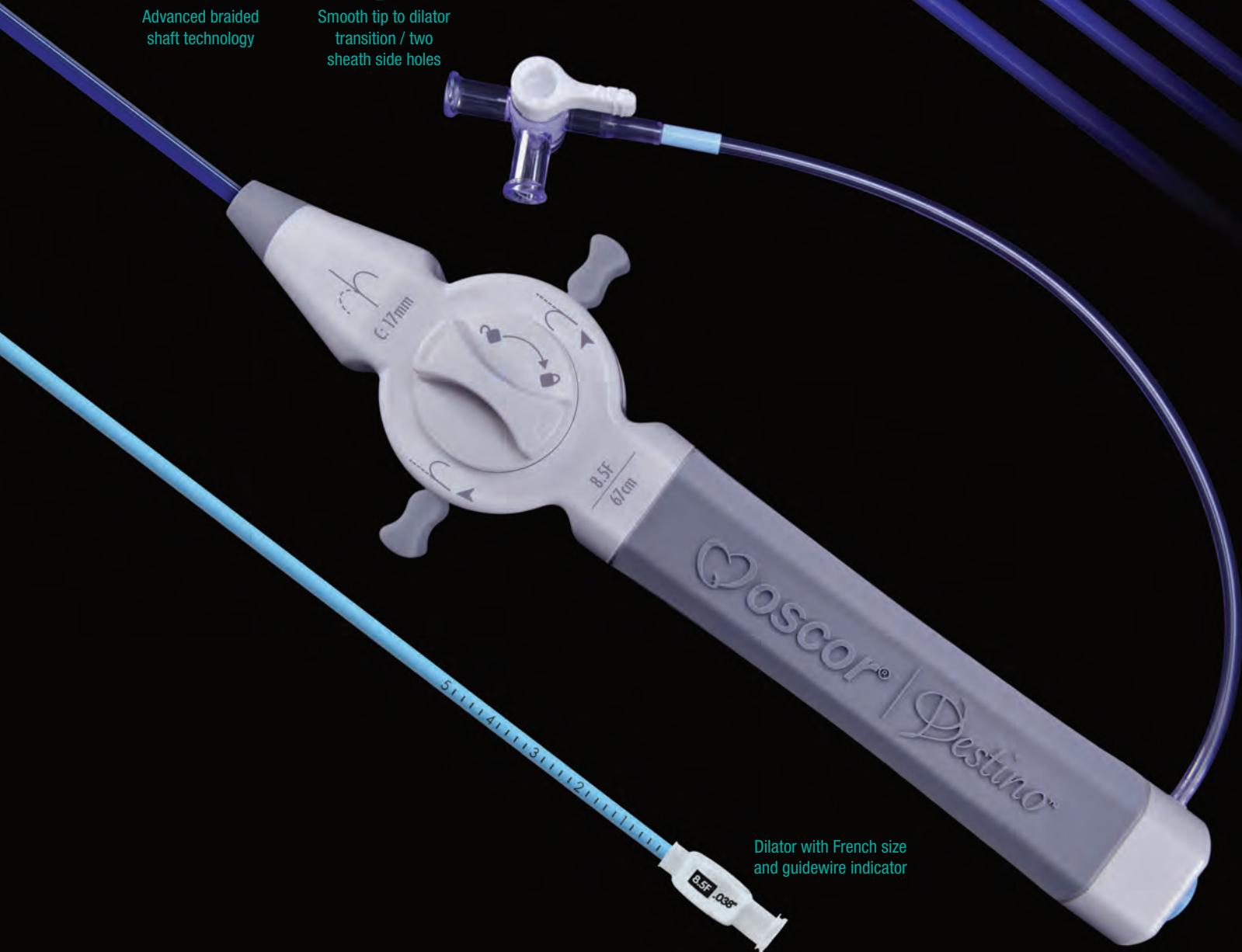
Destino®

Bi-directional Curve Options

50 mm



22 mm

17 mm



Dilator with French size and guidewire indicator

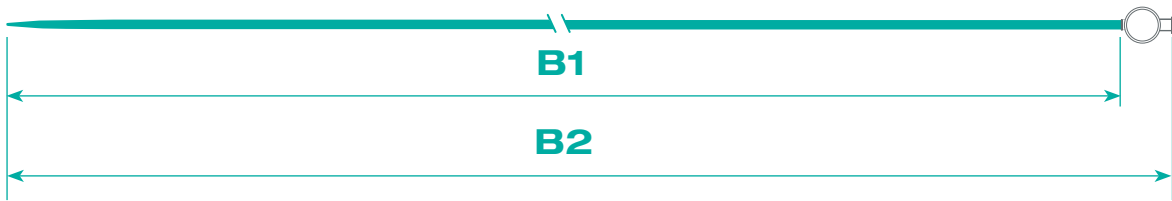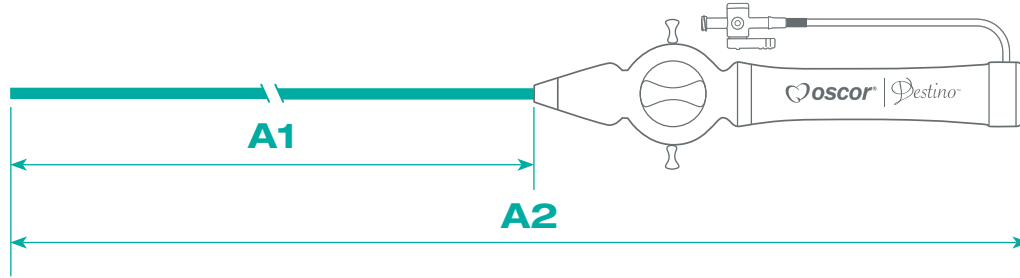Destino™ Finished Dimensions and Characteristics



Sheath Inner Diameter French Size	Tip Inner Diameter	Sheath Outer Diameter	A1 Usable Sheath Length	A2 Overall Sheath Length	B1 Usable Dilator Length	B2 Overall Dilator Length	C Compatible Transseptal Needle
8.5 F	0.112" (2.80 mm)	0.158" (4.10 mm)	71 cm	91 cm	94 cm	96 cm	98 cm
10 F	0.142" (3.60 mm)	0.185" (4.70 mm)	67 cm	87 cm	90 cm	92 cm	98 cm
12 F	0.168" (4.20 mm)	0.208" (5.30 mm)	67 cm	87 cm	90 cm	92 cm	98 cm

Guiding sheath supplied sterile. For single use only.

Standard lengths are noted above. Other lengths and curves available upon request.

Ordering Information

Model # Kit <i>(Straight Dilator)</i>	Model # Set <i>(Straight Dilator)</i>	Model # Kit <i>(50° Curved Dilator)</i>	Model # Set <i>(50° Curved Dilator)</i>	Sheath Inner Diameter <i>(French Size)</i>	Sheath Length	Sheath Curve Size
<i>8.5 F Series</i>						
DS857117GW	DS857117	DC857117GW	DC857117	● 8.5 F	71 cm	17 mm
DS857122GW	DS857122	DC857122GW	DC857122	● 8.5 F	71 cm	22 mm
DS857150GW	DS857150	DC857150GW	DC857150	● 8.5 F	71 cm	50 mm
<i>10 F Series</i>						
DS106722GW	DS106722	DC106722GW	DC106722	● 10 F	67 cm	22 mm
<i>12 F Series</i>						
DS126750GW	DS126750	DC126750GW	DC126750	● 12 F	67 cm	50 mm

Guiding sheath supplied sterile. For single use only.

Standard lengths are noted above. Other lengths and curves available upon request.

For additional information, e-mail us at sales@oscor.com

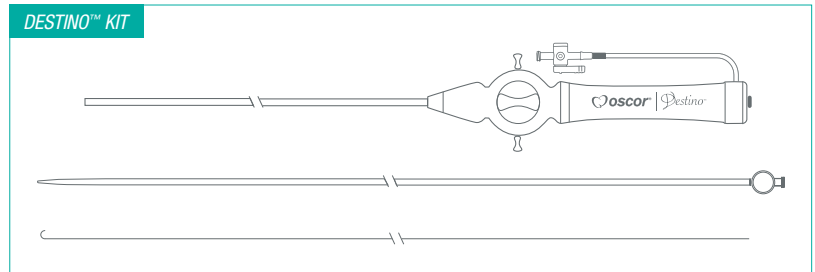
Products not available in all regions.

Destino™ Kit Contains:

- One (1) Destino™ steerable guiding sheath
- One (1) straight or 50° curved dilator
- One (1) guidewire OD: 0.035", L: 180 cm

Destino™ Set Contains:

- One (1) Destino™ steerable guiding sheath
- One (1) straight or 50° curved dilator



DS100-030121-00



OSCOR INC. | WORLD HEADQUARTERS

Toll Free: 1-800-726-7267
P: +1 727-937-2511

OSCOR EUROPE GMBH | GERMANY DIRECT SALES

Telefon kostenlos: 0800 00-67267
P: +49 (0) 211 586 786-00

E: sales@oscor.com
www.oscor.com



CAUTION: Federal Law restricts this device to sale by or on the order of a physician (U.S.A.). Please refer to the Instructions for Use for a complete listing of the indications, contraindications, precautions and warnings, where applicable.